Smart Neck Massager

User Manual

MODEL#4098/4353

Thank you for choosing this SKG product.
Before operating this unit, please read these instructions
completely and save them for future use.

Thank you for choosing our Smart Neck Massager. Please read this manual carefully before use in order to use the product correctly and take full advantage of its excellent performance. We believe it will bring convenience and comfort to your life.

Please keep this manual for your future reference. If any problems occur when using the product, please contact us. We will be glad to help you. Thank you!

Content

Content.....	1
Precautions.....	2
Technical parameters.....	5
Product diagram.....	6
Operation instruction.....	7
Tips.....	8
Daily maintenance.....	9
FAQ.....	9

This product is a household product. Don't use it for commercial purposes.

As we continually improve our products, we shall reserve the right to make technical changes without prior notice.

Precautions

Danger

- Please never use the product with the following medical electronic devices:
 1. Heart pacemakers and other implantable medical electronic devices;
 2. Artificial heart and lungs and other life-supporting medical electronic devices;
 3. Electrocardiographs and other medical electronic devices;
 4. Petrol stations, high dust workshops, gas and other inflammable and easy-explosive area are prohibition of use.

Warning

- If the patient applies high-frequency surgical device and the product on the body at the same time, the electrode of the product may cause burns and damage to the product.
- If the product is used near (1M) a short-wave or microwave therapy device, the output of the product may be unstable.
- The electrode should prevent contact with wounds and scars during treatment.
- When you use the product, don't place the electrode near the chest heart.
- The product is recommended for use by one user. When it's used by different users, clean the electrode on the product before using it.
- When you use the product, stick the electrode tightly and evenly to the skin.
- The electrode of this product has a heat-generating function. When the hot compress function is applied on the same part for a long time (more than 15 minutes), it may cause low temperature burns.
- If you feel unwell during or after use, stop using it and seek medical attention.
- This product contains a lithium battery. Don't disassemble, modify, flood the product, or throw it into fire, otherwise it may result in product damage and explosion.
- In the following cases, consult your doctor before use, otherwise it may cause accidents or cause physical discomfort.

1. Patients under medical treatment and feel discomfort.
 2. Patients with abnormal heart and brain nerves.
 3. Patients whose body temperature exceeds 38 degrees (fever period).
 4. Patients with infectious diseases.
 5. Patients with skin perceptual disorders or skin abnormalities.
 6. Patients under medication
 - 7, patients with abnormal blood pressure.
 8. Patients who are not free to express their wishes.
 9. Patients with sensory impairment due to high peripheral circulatory disorders such as diabetes.
 10. Patients with skin allergies to stainless steel electrodes.
- Usage under the following conditions is prohibited:
 1. Malignant tumor patients
 2. Patients with bleeding tendency
 3. Pregnant women, lying-in women
 4. Acute suppurative inflammation patients
 5. Heart pacemaker users
 6. Patients with metal implant
 7. Near the heart
 - Don't apply it near the heart, above the neck, head, mouth, genitals, or skin disorders and don't treat the soles of the feet simultaneously. In addition, don't compress the internal organs when using.
 - Don't use with other medical devices or more than 2 products.
 - Don't use for purposes other than the ones described in the instruction manual.
 - Don't use the electrode with a chipped or damaged surface.
 - Don't disassemble, repair or modify this product without authorization.

Note

- If you feel unwell due to an abnormality in this product, stop using it immediately.
- If your body or skin feel strange, consult your doctor and follow the doctor's advice.
- If you want to move the pad or place it on other person during massage, be sure to turn off the power before moving it.

- Don't treat while having electronic devices.
- If abnormal operation occurs, immediately switch off the power supply and remove it.
- Don't plug or unplug the USB charging cable with wet hands.
- Don't use it for children or people who have no expressiveness.
- Don't use it in a place with high humidity, such as a bathroom, or use it while bathing.
- Don't use it when you are sleeping.
- Don't stick the electrode to the affected area for a long time.
- When massaging, don't use the metal parts of belts, necklaces, etc., to touch the electrode.
- Don't apply the product on the same place for long time (more than 30 minutes).
- Please wear this product correctly and use it on the cervical spine.
- When unplugging the USB charging cable from the body, do not pull on the cable, but rather, pull the cable plug and pull it out.
- The product electrode has a surface heating function for users who are not sensitive to heat or subject to heat allergies.
- Please consult a doctor before use.
- When the product is working for a long time, the surface will be hot, which is normal.
- When the remote control battery is idle for a long time (more than 3 months), remove the battery.
- When the remote control battery is exhausted, replace it with a new battery of the same model as soon as possible. Battery model: CR2032.
- When replacing the remote control battery, do not install the pole reversibly.
- After the remote control battery is exhausted, dispose of it according to the relevant regulations in your place of residence.

Suggestion


- Please make sure that the host button is working properly.
- When you use it for the first time, or if the battery loses power after the product has been left idle for a long time, please perform the charging operation according to the instruction manual to ensure that it can work normally.
- When disposing of this product, please dispose of it in accordance with the relevant environmental protection regulations of the city, otherwise it may cause environmental pollution.
- It is recommended for people aged 16 and above.
- Don't use or store this product in the following locations:
 - Places exposed to direct sunlight
 - Places subject to high temperature and humidity
 - Places with water
 - Dusty places
 - Places near the fireworks
 - Places susceptible to vibration shocks
 - Strong electromagnetic places

Technical parameters


Model	4098	Rated voltage	5V 
Rated power	5W		
Standard	EN60335-1	EN60335-2-32	EN62233

Model	4353	Rated voltage	5V 
Rated power	5W		
Standard	EN60335-1	EN60335-2-32	EN62233


Product diagram

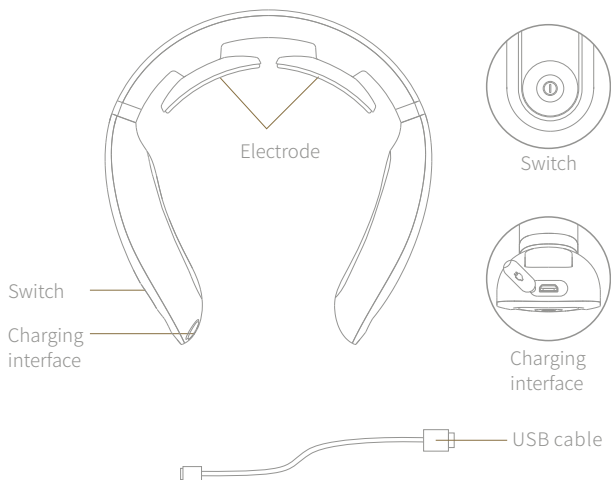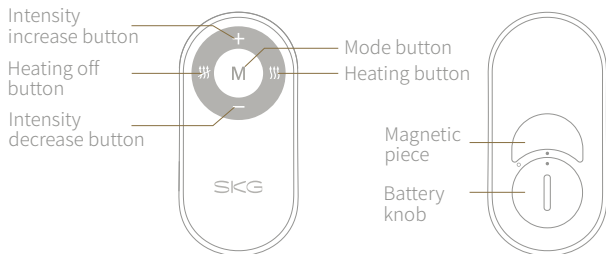Remote control

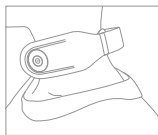Wearing method

Operation method

1. Before use, first wet the skin of the cervical spine with a wet towel, and then put the product around the neck, make the electrode close to the skin of the cervical spine, adjust the position of the product, turn on the host and adjust the appropriate massage intensity and mode.
2. Long press the host "  " button for two seconds, turn on the host (with the buzzer "beep" sound and prompt light), enter the default mode with the default level (0). Press the host button, add the intensity by one level with a "beep" sound. Until level 15, press the button and hear two "beep" sounds with no intensity change. Long press host "  " for two seconds when the machine is on, hear two "beep" sounds and turn off the products.
3. " **M** " on the remote control is the mode button. Press the mode button to switch the mode. There are three modes. After each press, it will send out "beep" sound and switch to next mode. The three modes can be switched in cycle (Note: daily mode, relax mode, exercise mode). Each time a mode is switched, the intensity returns to the default 0 level.
4. " **+** " on the remote control is the intensity increase button. After each press, it will send out "beep" sound and switch to higher level. The intensity is divided into 15 levels (0, 1, 2 ... 15). When the intensity is at level 15, when pressing " **+** ", it will send out two "beep" sound with same intensity at level 15.
5. " **-** " on the remote control is the intensity decrease button. After

each press, it will send out "beep" sound and switch to lower level. The intensity is divided into 15 levels (0, 1, 2 ... 15). When the intensity is at level 0, it will send out two "beep" sound with same intensity at level 15 after pressing " — ".

6. Long press the button " + " for two seconds to turn on the sound cues. Then press the button " — " for two seconds to turn off the sound cues.
7. "))) " on the remote control is the heating button. Press the heating button and start host heating function with a "beep" sound. The two electrode will be hot at the temperature of 37~43°C.
8. "))) " on the remote control is to heating off button. Press the heating off button and stop the heating function with two "beep" sounds.
9. The host is scheduled for 15 minutes. When the host is turned on, it will automatically shut down after 15 minutes.
10. When the host is turned on, the product will automatically shut down in 30 seconds when the product is not put on normally.

Product charging

1. When the host battery is running out, the power button lights up in orange. You can use the Micro-USB cable to plug into the other end of the host and plug in the charger with the Android head (the output is 5V1A/2A). When charging, the orange light flashes. When it is full, the blue light flashes.
2. When the remote control battery is replaced, a coin can be inserted into the battery knob slot. Turn it counterclockwise to the alignment position to open the replacement battery. Note that the battery positive pole is facing up. (Battery model: CR2032)

Tips

- The product will be automatically turned off after 15 minutes of operation. It is recommended to use it for no more than 30 minutes per day.
- The product should be stored in a dry and ventilated place after use.

- Do not place other weights on the product.
- Don't pull on the rubber plug at the USB charging port, as it may cause rupture and falling off.

Routine maintenance

- When cleaning and maintaining this product, please turn off the host, use a wet paper towel to clean the electrode. When it is particularly dirty, use a little medical alcohol (alcohol concentration 75%) to wipe and disinfect.
- Please wipe with neutral detergent and soft cloth during daily maintenance, do not use corrosive detergent and don't place it into water for soaking or washing.
- During the cleaning process, water or other liquids should be prevented from entering the body to avoid short circuit and electric shock during use.

FAQ

- If you have problems using the product, please refer to the following tips and solutions.
- If the product is still in a fault condition or cannot be used, turn off the power.
- Users should not disassemble or repair the product by themselves. Dismantling or repairing by themselves will cause damage to the product and cause the invalidity of warranty. The company will not be responsible for any consequences caused thereby.

Problems	Possible cause/treatment
Failure to start	<ul style="list-style-type: none"> ● The internal rechargeable battery is running low ● The remote control battery is dead. ● Please confirm whether the host product switch has been pressed.
Failure to charge or no charge light	Please check if the USB cable and host and phone power adapter are connected.

Problems	Possible cause/treatment
Normal display but no stimulation	<ul style="list-style-type: none"> ● Does the pad contact the skin directly? Adjust the position so that it is close to the skin. ● Is there hair or clothing, jewelry, etc. between the pads and the skin? Move away or remove items. ● Is the intensity level at 0? Increase the massage intensity.
Tingling sensation	<ul style="list-style-type: none"> ● Is the massager's electrode and neck skin closely contacted? adjust the position so that it' s stick to the skin. ● Is the massage time too long? It is recommended to massage for 15 minutes at a time, twice a day. ● Is the massage intensity too high? Adjust to low level ● Is the skin allergic to the electrode? Those with severe allergies need to be desensitized first. ● Is the climate cold and the skin dry? Use a wet towel to moisten the skin and use it. ● Does the skin have a wound? Do not put the electrodes on the wound for massage.
Weak simulation	<ul style="list-style-type: none"> ● Is the massager's electrode and neck skin closely contacted? adjust the position so that it' s stick to the skin and adjust intensity. ● Is the electrode dirty or oily? Clean the electrode surface. ● Is the battery running low? Please charge if the battery is low. ● Is massage intensity level too low? Increase massage intensity.

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